



ID	Compound	Rt (min)	ID	Compound	Rt (min)	ID	Compound	Rt (min)	ID	Compound	Rt (min)
1	PCB4 / PCB10	25,10	13	PCB31	35,01	25	PCB40	40,85	38	PCB85	46,58
2	PCB7 / PCB9	26,96	14	PCB28	35,12	26	PCB67 / PCB100	41,28	39	PCB136	46,79
3	PCB6	27,80	15	PCB20 / PCB33 / PCB53	35,85	27	PCB63	41,69	40	PCB110	46,98
4	PCB5 / PCB8	28,25	16	PCB22 / PCB51	36,45	28	PCB74	41,98	41	PCB151 / PCB82	47,88
5	PCB19	29,84	17	PCB45	36,84	29	PCB70	42,23	42	PCB135	48,19
6	PCB30 (IS)	30,48	18	PCB46	37,44	30	PCB66 / PCB95	42,57	43	PCB107	48,51
7	PCB18	31,55	19	PCB52 / PCB69	37,73	32	PCB56 / PCB60	43,79	44	PCB149 / PCB123	48,74
8	PCB15 / PCB17	31,71	20	PCB49	38,08	33	PCB101 / PCB90	44,26	45	PCB118	48,87
9	PCB24 / PCB27	32,36	21	PCB47 / PCB48 / PCB75	38,39	34	PCB99	44,67	46	PCB134	49,52
10	PCB16 / PCB32	33,02	22	PCB44	39,30	35	PCB83	45,50	47	PCB146	50,05
11	PCB26	34,37	23	PCB37 / PCB42 / PCB59	39,58	36	PCB97	45,86	48	PCB153	50,50
12	PCB25	34,58	24	PCB41 / PCB64	40,25	37	PCB87 / PCB115	46,24	49	PCB132	50,71

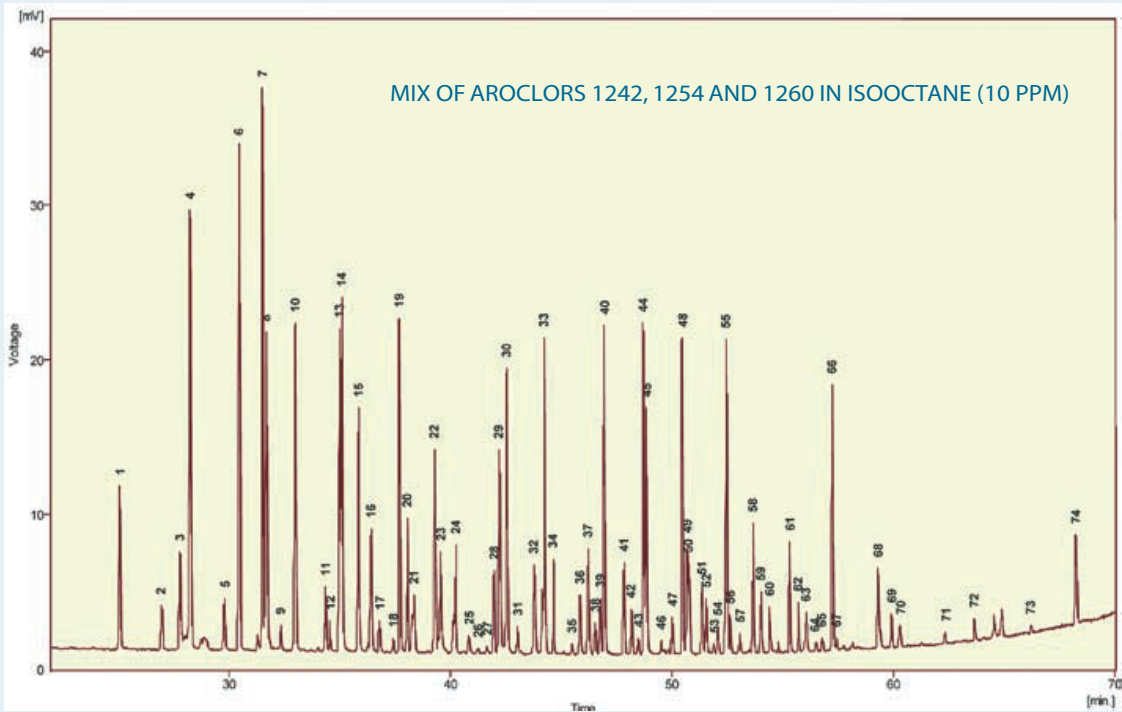

50	PCB105	50,79
51	PCB141	51,40
52	PCB179	51,57
53	PCB130	51,93
54	PCB137 / PCB176	52,08
55	PCB138 / PCB160 / PCB163	52,48
56	PCB158	52,67
57	PCB129 / PCB178	53,11
58	PCB187	53,69
59	PCB183	54,05
60	PCB128	54,43
61	PCB174	55,35
62	PCB177	55,73
63	PCB171 / PCB156	56,07
64	PCB201 / PCB157 / PCB173	56,51
65	PCB172	56,80
66	PCB180	57,27
67	PCB193	57,47
68	PCB170 / PCB190	59,33
69	PCB199	59,94
70	PCB196 / PCB203	60,34
71	PCB195 / PCB208	62,38
72	PCB194	63,68
73	PCB206	66,23
74	PCB209 (IS)	68,25

cronus

AROCLOR 1242 IN ISOCTANE (10 PPM)

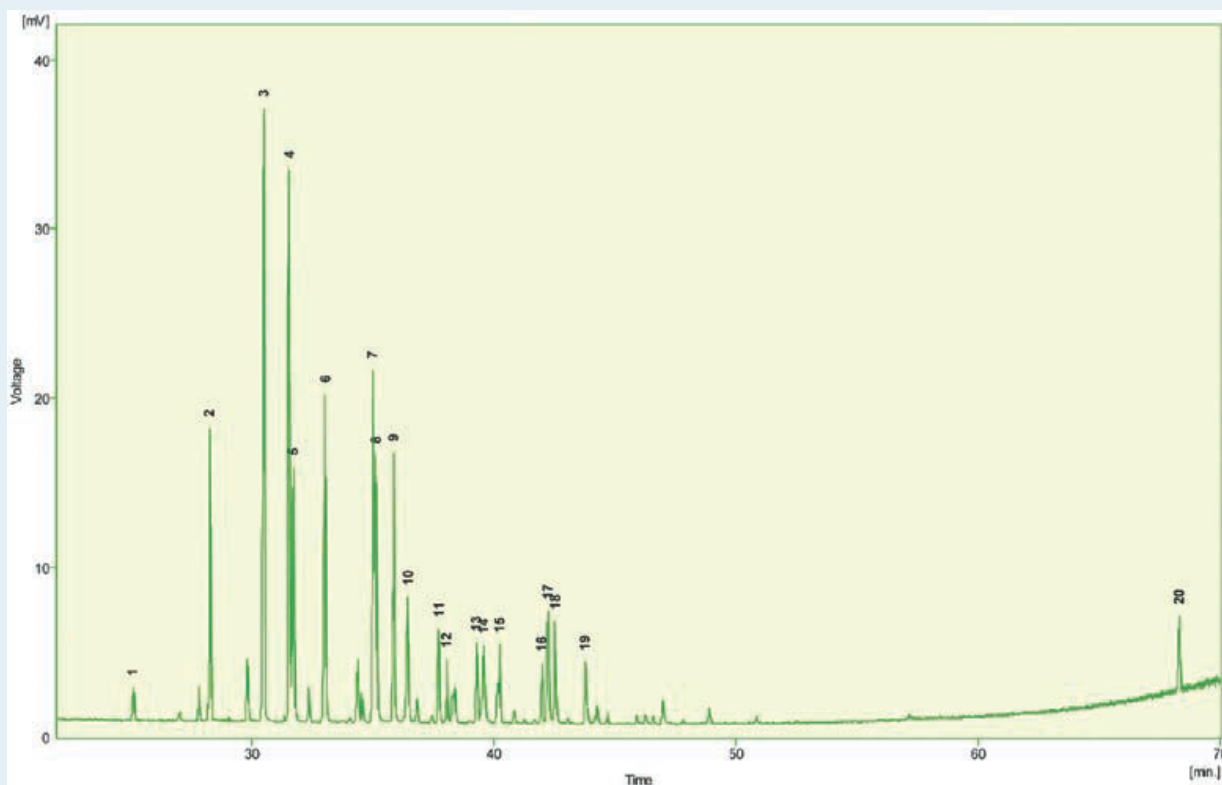

Peak Name	ID	Compound
1	PCB4	PCB4
2	PCB8	PCB8
3	PCB30 (IS)	PCB30 (IS)
4	PCB18	PCB18
5	PCB15/PCB17	PCB15/PCB17
6	PCB16/PCB32	PCB16/PCB32
7	PCB31	PCB31
8	PCB28	PCB28
9	PCB33	PCB33
10	PCB22	PCB22
11	PCB52	PCB52
12	PCB49	PCB49
13	PCB44	PCB44
14	PCB37/PCB42/PCB59	PCB37/PCB42/PCB59
15	PCB41/PCB64	PCB41/PCB64
16	PCB74	PCB74
17	PCB70	PCB70
18	PCB66	PCB66
19	PCB56/PCB60	PCB56/PCB60
20	PCB209 (IS)	PCB209 (IS)

cronus