

SEPARATION OF SUGARS (AS ALDITOL ACETATES)

Column: **CGB-1701**, P/N CG-130212


Dimensions: 15m x 0.25mm x 0.25 μ m

Injection: 1 μ L Sugars standard, split

Carrier gas: H₂, 6 psi (41.3 KPa)

Oven temperature: 180°C @ 4°C/min to 215°C

Detector: FID, 250°C

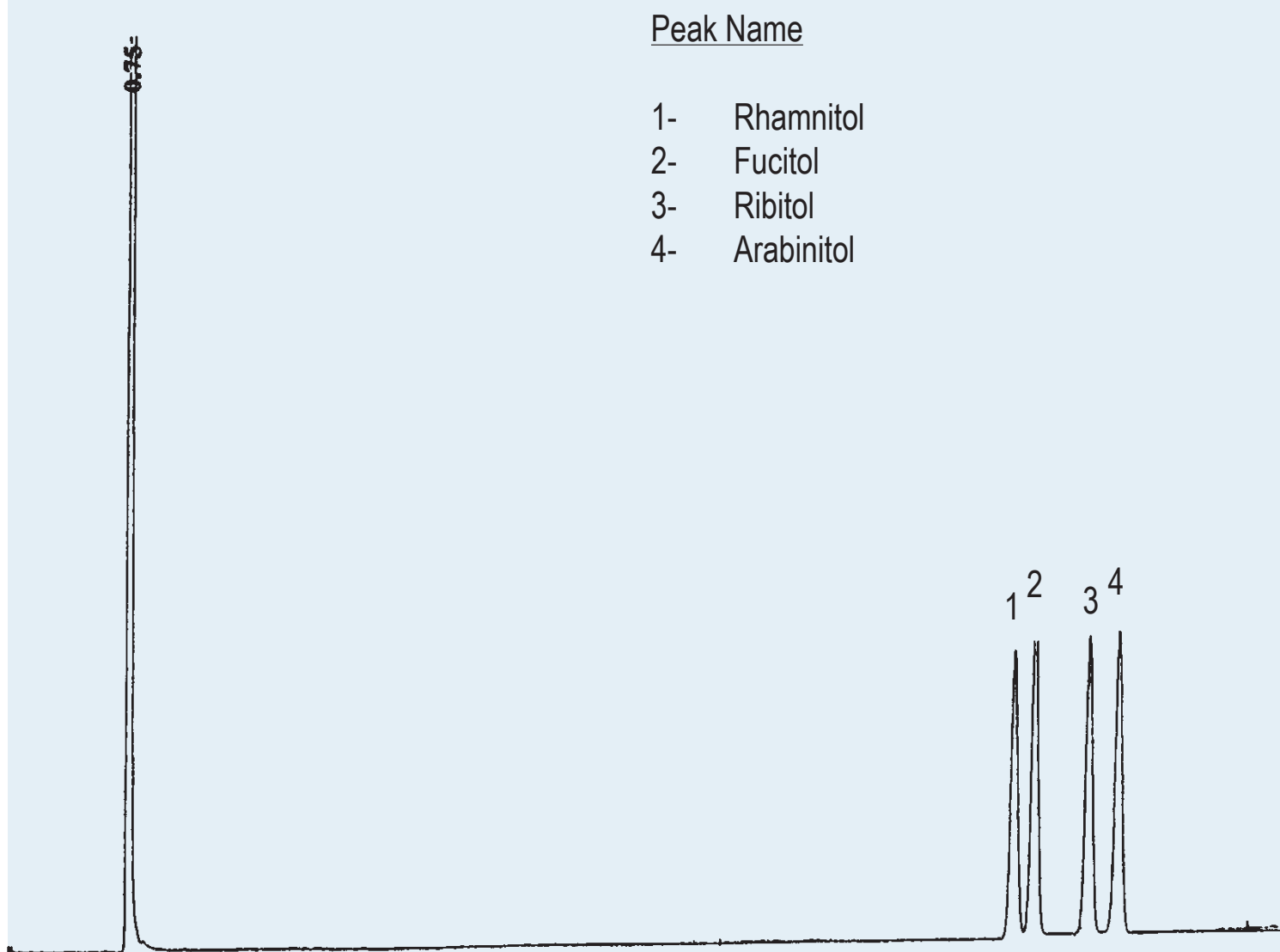

Peak Name

- 1- Rhamnitol
- 2- Fucitol
- 3- Ribitol
- 4- Arabinitol