
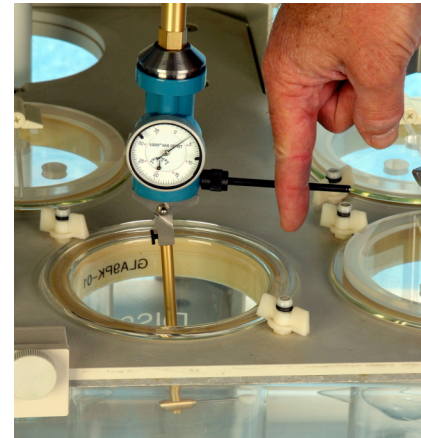


# CNTGAG-MI

## Multi-Instrument, Alignment Gage

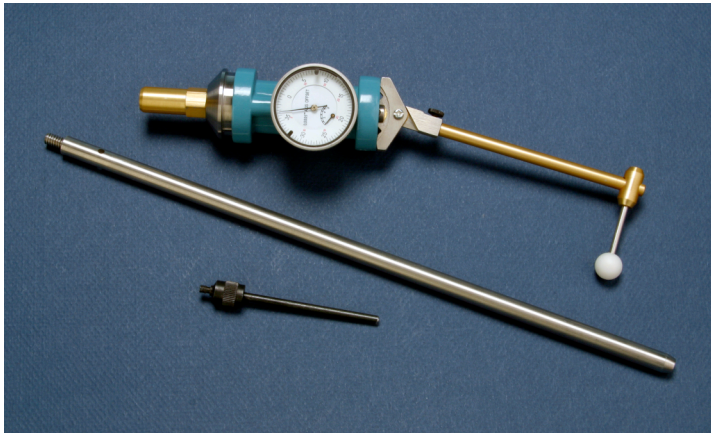

USP 711 states: “the shaft is positioned so it’s axis is no more than 2 mm at any point from the vertical axis of the vessel.” Prior art instruments with dial indicators which rotate are difficult if not impossible to read. CNTGAG-MI has a large dial indicator which does not rotate and is easy to read at all vessel positions.

- ✓ *Quickly and easily mounts to all vessel positions*
- ✓ *Works with all instruments which have access to the vessels in the operating position*
- ✓ *Non rotating, easy to read Dial Indicator*
- ✓ *High resolution – One division = .05 mm*
- ✓ *Range : More than 4 mm*
- ✓ *NIST Traceable Calibration*

### Ordering Information:

Catalog Number	Description
CNTGAG-MI	Alignment Gage, Multi-Instrument